

INSTRUCTIONS FOR

LUD Press & Holder for BD Tamper-Tuf™ Vials with Caps (#18737 & 18738)



Step 1. Place the vials into the holes of the holder.



Step 2. Continue placing the vials in the holder so that all of the openings are aligned.



Step 3. Start at one end of the holder and fill the first vial. Continue moving from one vial to another until all of them are filled.



Step 4. Slide the holder into the channel base of the LUD Press. Make sure that the tabs are facing away from the press.



Step 5. Using your finger, push the tab down on the vial about to be sealed.



Step 6. Position the tab under the press of the and push down until you hear the tab snap closed.

Call Free 1.800.848.1633 • Fax Free 1.800.447.2923

Web: GoHCL.com • Email: hcl@GoHCL.com

INSTRUCTIONS FOR

LUD Press & Holder for BD Tamper-Tuf™ Vials with Caps (#18737 & 18738)



Step 7. Slide the holder forward to the next vial.

Call Free 1.800.848.1633 • Fax Free 1.800.447.2923
Web: GoHCL.com • Email: hcl@GoHCL.com

INSTRUCTIONS FOR

LUD Press & Holder for Easy Fill Vials with Plugs (#18737 & 18738)



Step 1. Remove the top plate from the holder then place the vials top down onto the holder.



Step 2. Place the top plate over top of the vials locking them into place for the filling process.



Step 3. Start at one end of the holder and fill the first vial. Continue moving from one vial to another until all of them are filled.



Step 4. Remove the fitting from the base of the LUD Press. Place it in the barrel of the press and hand tighten it into place.



Step 5. Slide the holder into the channel base of the LUD Press.



Step 6. Position the plugs into the end of each vial.

Call Free 1.800.848.1633 • Fax Free 1.800.447.2923
Web: GoHCL.com • Email: hcl@GoHCL.com

INSTRUCTIONS FOR


LUD Press & Holder for Easy Fill Vials with Plugs (#18737 & 18738)



Step 7. Position the plug under the press and push down until you feel the plug push into place.



Step 8. After all the vials are sealed remove the top plate.

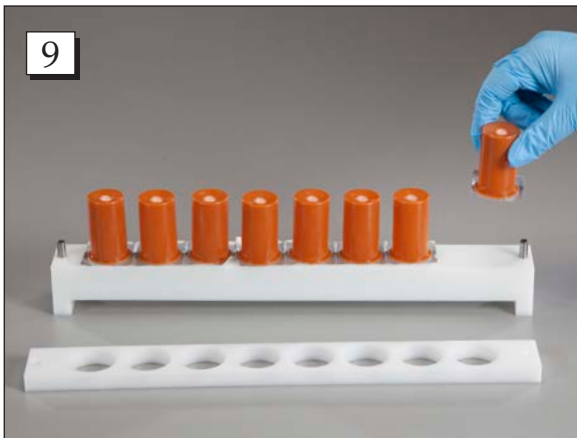

Step 9. Remove the sealed vials for use.

Call Free 1.800.848.1633 • Fax Free 1.800.447.2923
Web: GoHCL.com • Email: hcl@GoHCL.com