


ELECTRONIC BALANCE CALIBRATION MASSES

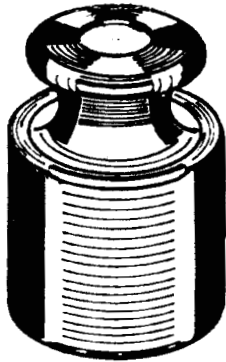

Directions for use
and Maintenance

Please read this
brochure before
using (or even
touching) your
Calibration
Masses

CALIBRATION IS
easier and more accurate.

Your new CALIBRATION MASS has been designed for highly precise calibration and adjustment of Electronic Balances.

The self-hinged carrying case provides a compartment to protect the mass when not in use.

USE

Electronic balances are normally calibrated by placing a known mass on the platform and adjusting the displayed weight to match that of the calibration mass. Refer to the Balance Instruction Manual for specific procedures.

CARE AND MAINTENANCE

Your Electronic Balance calibration mass(es) must be handled carefully and the following precautions are mandatory:

- **PROTECT** against **DAMAGE**
- **DO NOT TOUCH** except when using handling glove(s). Oil from hands will alter the tolerance of the mass.
- **KEEP CLEAN** and **DRY**. Any build up of foreign matter will alter the tolerance of the mass.
- **DO NOT DISTURB** the sealing plug or knob.

Call Free: 1 800 848 1633

Fax Free: 1 800 447 2923

Email: hcl@hcl-intl.com